

CR-750 CR-750S



CR-750

Operable via foot controller



CR-750S

Controllable from the ST-50 Smart Table

Designed to support clinical workflow

Serving Your Vision

TAKAGI ophthalmic devices, proudly made in Japan since 1955

This motorised patient chair combines patient comfort with ease of operation for medical staff. Smooth height adjustment enables easy posture changes during examinations, while a design that blends with the clinic environment supports a cohesive clinical setting.

- 150 mm vertical stroke accommodates a wide range of patients
- Wide legroom and smooth operation
- Designed to meet strict RoHS 2 standards



learn more

Patient comfort and efficient examinations


The 150 mm vertical stroke supports comfortable positioning for patients of various body types. In addition, the optional sliding mechanism allows patients to be positioned closer to examination equipment without moving the chair, helping to streamline examinations.



Staff-friendly operation

In addition to operation via the UP/DOWN switch on the back, the CR-750 can be operated using the standard foot controller, and the CR-750S can be operated from the ST-50 Smart Table. This allows for smooth and precise adjustments during use, reducing the workload for medical staff.

**The CR-750 cannot be connected to the ST-50.*

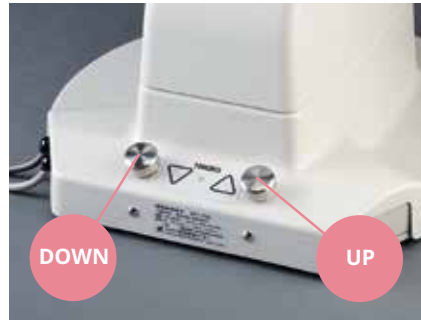

Complies with RoHS 2

The CR-750 and CR-750S are compliant with the EU RoHS 2 Directive (2011/65/EU) and free from restricted hazardous substances.

**The CR-750S is compatible with the Smart Table (ST-50).*




Specifications

Model name	CR-750 Motorized Patient Chair CR-750S Motorized Patient Chair
Elevation stroke	150 mm
Load capacity	up to 150 kg
Input power	AC 100-230 V, 50-60 Hz
Power consumption	120 VA
Weight	26.5 kg (CR-750), 25 kg (CR-750S)

Dimensions

Dimensions (mm)

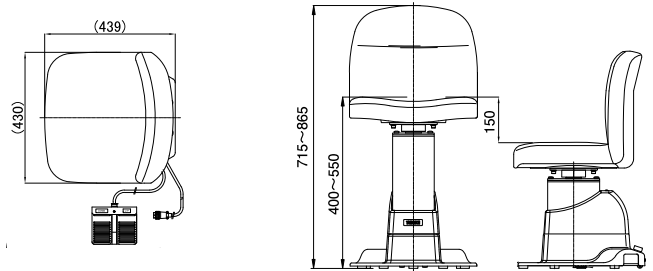

**The dimensions of CR-750 and CR-750S are identical.
CR-750S does not include the foot controller.

Description and appearance as detailed in this brochure may be subject to change as improvements are made to products. Colours of product(s) in brochure may slightly differ from actual product(s) due to lighting situation when photographed or print conditions.

To ensure safe and correct handling please read the user manual before using

Takagi Ophthalmic Instruments Europe Ltd.

Citylabs 1.0
(Former Royal Eye Hospital Building)
Nelson Street
Manchester, M13 9NQ
United Kingdom
T: +44 (0)161 273 6330
W: <https://www.takagieurope.com>



Takagi Seiko Co., Ltd.



takagieurope.com

Your local distributor