

ST-26T/ST-26H



Wheelchair friendly motorised table

Stylish and versatile medical grade electric instrument table, considering extreme safety with quiet operation for the most demanding medical environments.



- Versatile Height Accommodation
- Space-Efficient Footprint
- High Load Capacity



Wheelchair friendly

The large underfoot space of the ST-26 allows a patient in a wheelchair to be examined quickly and comfortably.

The ST-26T with a T-shaped base has its column at the side of the table. The ST-26H with a H-shaped base has its column at the center of the table.

* Wheelchairs with a front wheel outer width of 550mm or less are adaptable.

* It is recommended to fold up the footsteps of the wheelchair during examination.



Careful consideration of cable positioning

Access hole is placed on the table top to pass a power cable through, and three auxiliary power outlets under the table top allow to power a full range of products without any need to run extra cables to a standard mains power supply.

The mains power cable for the ST-26H extends from the side of the base eliminating interference with patient's feet.



Heavy load capability

With a load capacity of 80 kg, the ST-26H is well designed to withstand heavy equipment such as MT-325UD projection perimeter.



Large elevation stroke for all patient heights

Minimum height of the table is now approximately 100 mm lower than that of the previous model meaning small children can be easily examined. Without changing maximum height from the previous models, both the ST-26T and ST-26H have a 300 mm elevation stroke to accommodate small to large patients sufficiently.

Space saving design

While the size of the table top for the ST-26T remains the same as our previous model, the design change of the elevation unit achieved an overall width reduction of 14%, enabling ST-26T to be used efficiently in a limited space.

Various kinds of table tops

Various shapes of table tops for the ST-26H are available to choose according to user requirements and applications.


* The optional 600mm x 400mm table top is recommended.

* Table top for the MT-325UD Perimeter is included as standard with MT-325UD.

* Please contact our sales department for full details on available table tops.

Dimensions

Dimensions (mm)

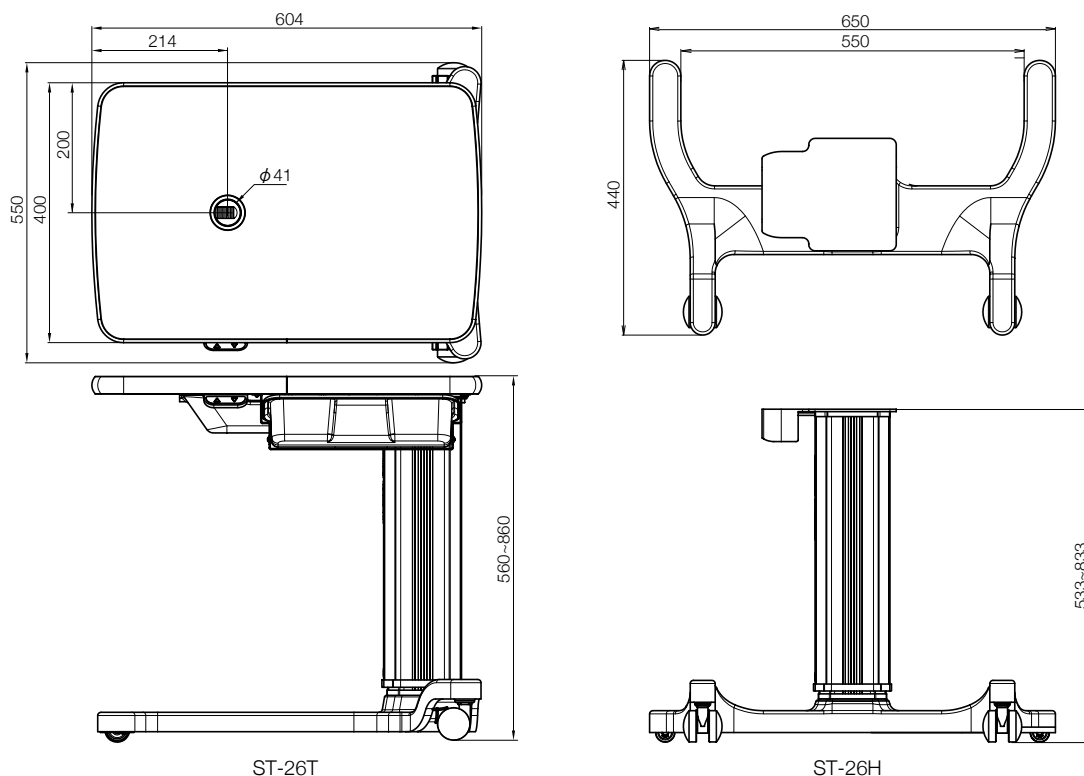

*The table top for ST-26H is optional.

Specifications

Model name	ST-26T Electric Instrument Table	ST-26H Electric Instrument Table
Minimum height	560 mm	533 mm *Not including table top
Maximum height	860 mm	833 mm *Not including table top
Elevation stroke	300 mm	
Allowable load	30 kg	80 kg
Power	AC 100-240 V, 50/60 Hz	
Power consumption	950-2150 VA	
Auxiliary electrical outlets	3 outlets	
Weight	20.2 kg	15 kg
Size of base	604 (W) × 550 (D) mm	650 (W) × 440 (D) mm
Size of table top	600 (W) × 400 (D) mm	-

Description and appearance as detailed in this brochure may be subject to change as improvements are made to products.

Colours of product(s) in brochure may slightly differ from actual product(s) due to lighting situation when photographed or print conditions.

 To ensure safe and correct handling please read the user manual before using

Takagi Ophthalmic Instruments Europe Ltd.

Citylabs 1.0
(Former Royal Eye Hospital Building)
Nelson Street
Manchester, M13 9NQ
United Kingdom
T: +44 (0)161 273 6330
W: <https://www.takagieurope.com>



Takagi Seiko Co., Ltd.



takagieurope.com

Your local distributor