

ST-50



Smart table with universal design

A smart table designed for easy slit lamp use.
Accommodates patient movement and facilitates ease of use for the clinician, with a compact design that fits neatly into the examination room.

- Comfortable for both doctors and patients
- Flexible operation supports smooth examinations
- Space-saving smart design



learn more

Design optimised for clinical environments

Designed with consideration for the examination environment of both doctors and patients, contributing to a comfortable clinical setting.




Safety switch

The underside of the table is fully covered to prevent exposure of mechanical components. The entire cover functions as a safety switch, helping to prevent pinching and accidental contact. It also helps prevent clothing from being stained by machine oil and supports cleanliness in the examination room.

Thin tabletop design

The thin design ensures ample space beneath the table, helping reduce the risk of contact with hands or thighs and contributing to reduced physical strain during examinations.

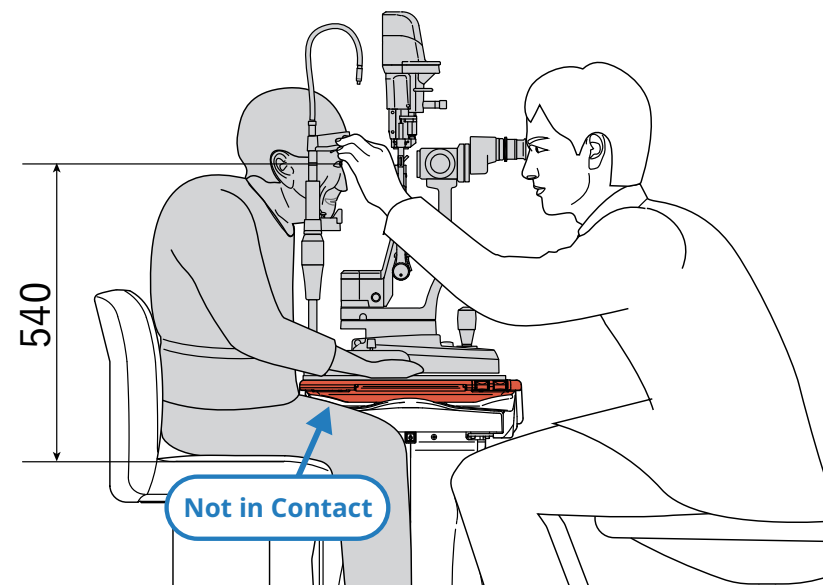

Easy operation for smooth alignment of the slit lamp and patient

Flexible operation and a design suited to various examination scenarios help improve operability during examinations.



Chair elevation switch

The table provides power and vertical movement control for the patient chair, supporting a smooth clinical workflow.

**Compatible model: CR-750S Motorized Patient Chair and UN-21 Full Reclining Chair.*

**CR-750S comes standard with a cable. Please see the Accessories on the back page for UN-21 cable.*



Elevation stroke

With a minimum height of 620 mm and a maximum height of 920 mm, the table accommodates examinations for patients with a wide range of body types.

Auto slide / manual slide

Automatic or manual slide methods can be selected depending on the examination situation and patient condition.

A smart design that complements clinic interiors

Designed to facilitate easy placement of peripheral equipment, enabling efficient layout configuration.

Slit lamp installation kits

Supports installation of slit lamps from various manufacturers, allowing flexible setup options.

**Please contact our sales department for full details.*

Easy-to-organise design

The sliding tabletop and cable outlets support efficient arrangement of peripheral equipment and help maintain an organised clinical space.



Combination with UN-21 Full Reclining Chair

In medical environments requiring large and/or reclining patient chairs, the UN-21 Full Reclining Chair can be connected to the table.

**Please see the Accessories on the back page for UN-21 cable.*




Specifications

Model name	ST-50 Smart Table
Minimum height	620 mm
Maximum height	920 mm
Slide stroke	500 mm
Allowable load	30 kg
Power	AC 100-230 V, 50-60 Hz
Power consumption	1500 VA
Weight	103 kg (including balancing weights exclusive of options)

Dimensions

Dimensions (mm)

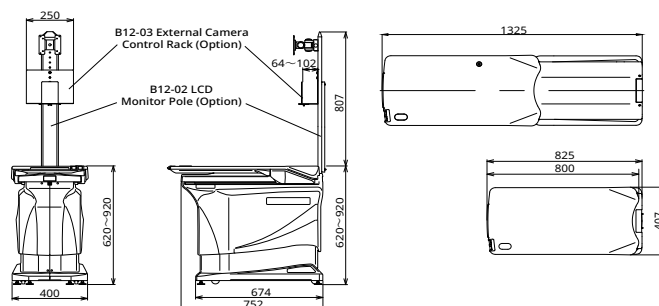

Accessories

Flexible customisation with a rich variety of options

1	B12-02	LCD Monitor Pole	An LCD monitor can be attached
2	B12-03	External Camera Control Rack	A camera control unit can be stored
3	B12-14	Scope Holder type B	For left-handed
4	B12-15	Scope Holder type A	For right-handed
	U23-01	Connection Cable	UN-21 can be connected

*Outer and inner covers are included with the scope holder.


*The photo on the right is an example of a system desk.



Layout

Available in two types: Type A and Type B.

These types are factory options and cannot be changed after delivery. Please specify the appropriate type according to your examination room layout when ordering.

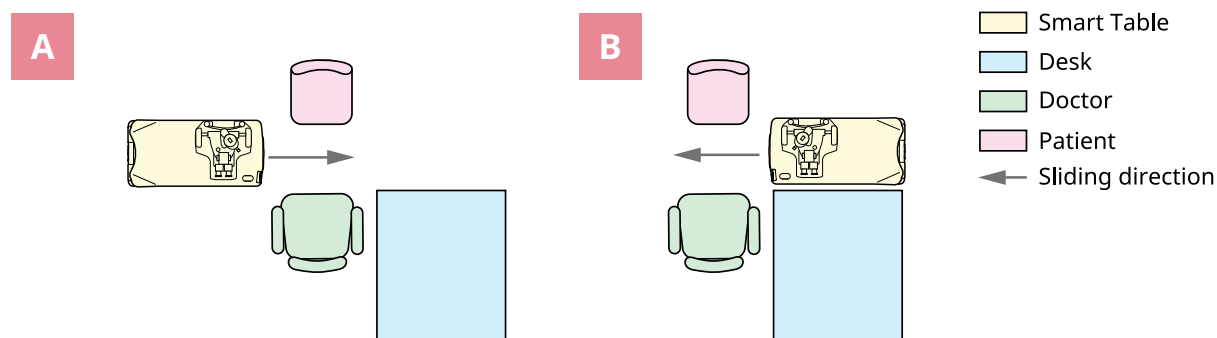

Description and appearance as detailed in this brochure may be subject to change as improvements are made to products.

Colours of product(s) in brochure may slightly differ from actual product(s) due to lighting situation when photographed or print conditions.

 To ensure safe and correct handling please read the user manual before using

Takagi Ophthalmic Instruments Europe Ltd.

Citylabs 1.0
(Former Royal Eye Hospital Building)
Nelson Street
Manchester, M13 9NQ
United Kingdom
T: +44 (0)161 273 6330
W: <https://www.takagieurope.com>



Takagi Seiko Co., Ltd.



[takagieurope.com](https://www.takagieurope.com)

Your local distributor