

UN-15



Hydraulic up/down chair with pneumatic assisted full recline capability

Built for daily use, designed for patient care

A pneumatic-assisted manual recline mechanism uses patient's weight for easy reclining adjustments.

Recline capability includes a comfortable flat position, or any desirable position in between.

Single-handed headrest operation provides smooth patient headrest adjustment, and a unique lock indicator ensures that the headrest is locked firmly and securely in position.

TAKAGI standard fail-safe sealed hydraulic system ensures smooth, long-term operation.

Rear-positioned, foot-activated rotation lock is easily accessible from both sides.



Single-handed headrest




Foot-activated rotation lock

Specifications

Model name	UN-15 Ophthalmic Chair
Reclining angle	107.5° to 172.5°
Rotating angle	255°
Reclining	By power hydraulic pump motor
Allowable load	150 kg
Power supply	AC 100, 120, 230 V, 50/60 Hz
Power consumption	800 VA
Weight	110 kg

Dimensions

Dimensions (mm)

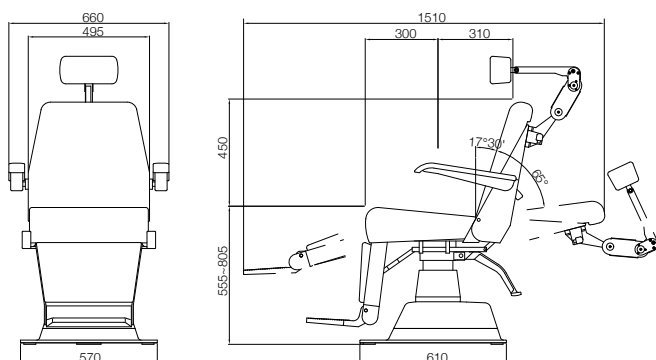

Combined with the SN-500 Ophthalmic Stand

Description and appearance as detailed in this brochure may be subject to change as improvements are made to products.

Colours of product(s) in brochure may slightly differ from actual product(s) due to lighting situation when photographed or print conditions.

 To ensure safe and correct handling please read the user manual before using.

Takagi Ophthalmic Instruments Europe Ltd.

Citylabs 1.0
(Former Royal Eye Hospital Building)
Nelson Street
Manchester, M13 9NQ
United Kingdom
T: +44 (0)161 273 6330
W: <https://www.takagieurope.com>



Takagi Seiko Co., Ltd.



learn more



[takagieurope.com](https://www.takagieurope.com)

Your local distributor