

UN-21



Full reclining patient chair with memory mechanism
and a foot controller

In combination with an operating microscope, the chair can be served as a bed for examination and treatment.

(An optional Connecting Arm (U21-02) is required to install an operating microscope.)

An AUTO-RETURN switch automatically returns the chair to its original position.

A standard memory mechanism can preselect one desirable tilt angle.

A standard foot switch controls lifting/lowering, back/forward tilting, and stopping positions of the chair.

The environmentally friendly leather is used for the compliance with RoHS Directive.

Leather fabric physical properties :


Alcohol resistance: Soaked in 70% ethanol water solution for 24 hours. Hypochlorous acid resistance: Soaked in 0.5% hypochlorous acid for 24 hours. Flameproof: Conforms to the flameproof product performance test standards of the Japan Fire Retardant Association. Antibacterial, light resistance, and durability conform to JIS (Japanese Industrial Standards).

Specifications

Model name	UN-21 Full Reclining Chair
Rotation	30° clockwise, 225° counterclockwise from step side
Elevation	By power hydraulic pump motor
Elevation stroke	250 mm
Reclining	By power hydraulic pump motor
Reclining angle	65° from upright position
Allowable load	150 kg
Power supply	AC 120, 230 V, 50/60 Hz
Power consumption	800 VA
Weight	121 kg

Dimensions

Dimensions (mm)

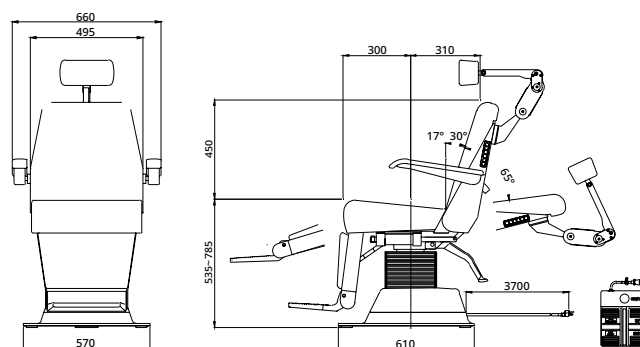

Description and appearance as detailed in this brochure may be subject to change as improvements are made to products.

Colours of product(s) in brochure may slightly differ from actual product(s) due to lighting situation when photographed or print conditions.

To ensure safe and correct handling please read the user manual before using.

Takagi Ophthalmic Instruments Europe Ltd.

Citylabs 1.0
(Former Royal Eye Hospital Building)
Nelson Street
Manchester, M13 9NQ
United Kingdom
T: +44 (0)161 273 6330
W: <https://www.takagieurope.com>



Takagi Seiko Co., Ltd.



learn more



takagieurope.com

Your local distributor

Combination with ST-50



When connected, the ST-50 provides power and vertical movement control for the UN-21.

Please see the ST-50 brochure for details.